


## Submittal Data Sheet

### Features


- Adjustable line regulator 0 – 125 PSIG or 0 – 200 PSIG (Nitrogen)
- Flexible stainless braided pigtails CGA gas specific pigtail connections are provided standard for all gas services, Copper pigtails are available as an option
- Complete system cleaned for use with Oxygen
- Cylinder contents gauge and delivery pressure gauge included
- Outlet union provided – 1/2 Male NPT
- Wall mounting bracket included
- Maximum inlet pressure 3,000 PSIG
- Available with pressure switch port 1/2 female NPT
- Optional pressure switch set to 150 PSI
- Optional remote audio visual alarm preprogrammed to “OXYGEN LOW SWITCH MANIFOLD” signal

### Specifications

Ideal for veterinary applications where manual switching of cylinders and a possible interruption of gas service are allowable. Systems are available for Nitrous Oxide, Carbon Dioxide, Nitrogen, Medical Air and Oxygen.

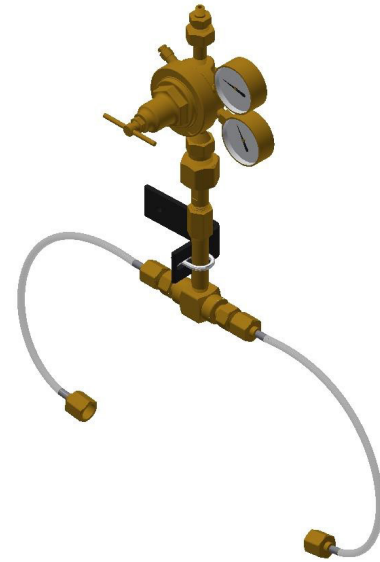


Model PX-TSD-540 with flexible pigtails shown above

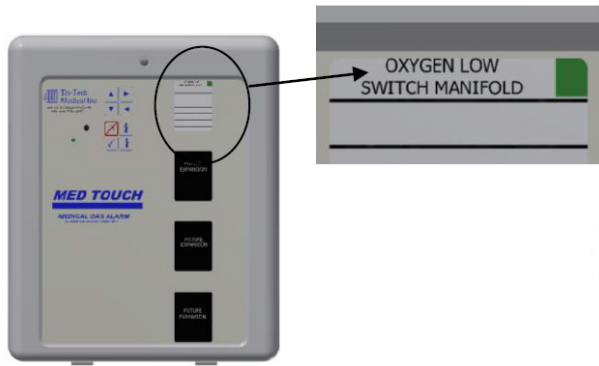

Model PX-ST2UPFFF-02 shown above  
Optional remote audio visual 8 point alarm must be used with optional pressure switch



Model PX-SDWPHP-580 with PX-PS-100-1000 optional pressure switch shown above

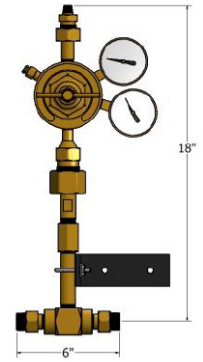

Dimensional information  
Allow 20" of wall space for two gas cylinders.


## Ordering Information

Easy to use modular ordering system.

Fill in the blanks to specify the TSD that meets your needs.

PX- 

TSD	P	D	GGG
-----	---	---	-----

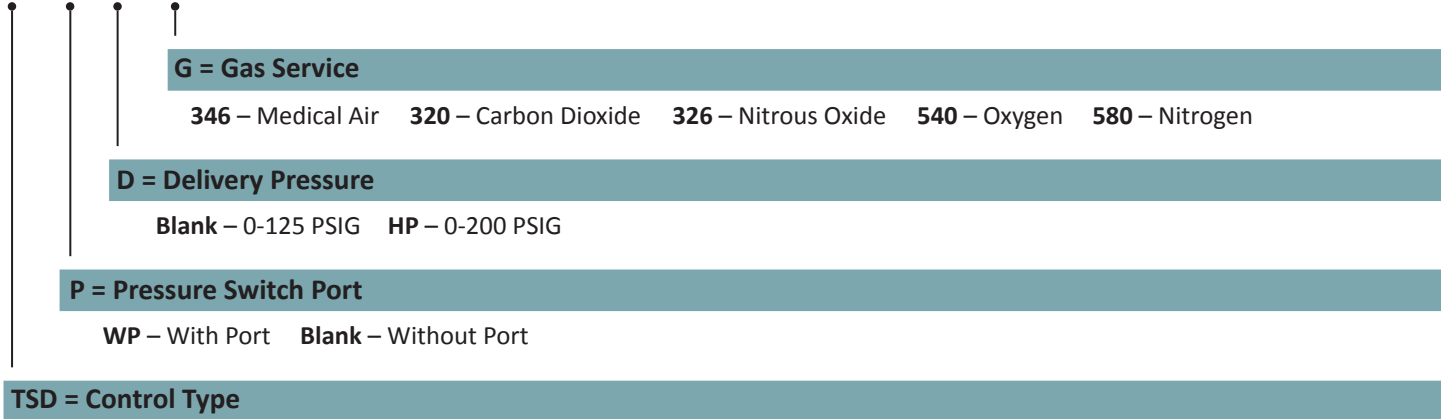

TSD – Type of Manifold

**Examples:** PX-TSD-540 = TSD system without pressure switch port, 0 – 125 psig, for Oxygen  
PX-TSDWPHP-580 = TSD system with pressure switch port, 0 – 200 psig, for Nitrogen