
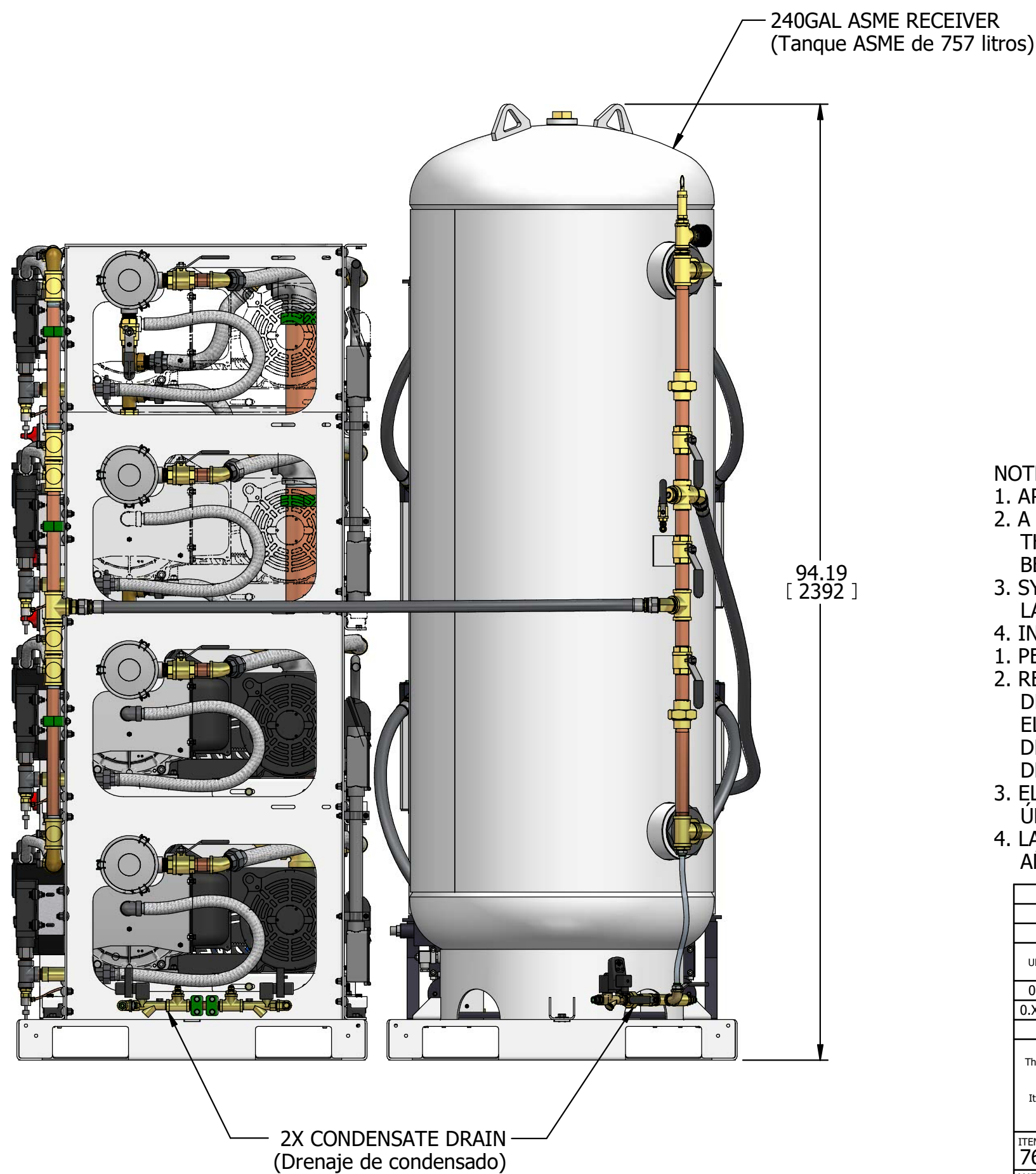
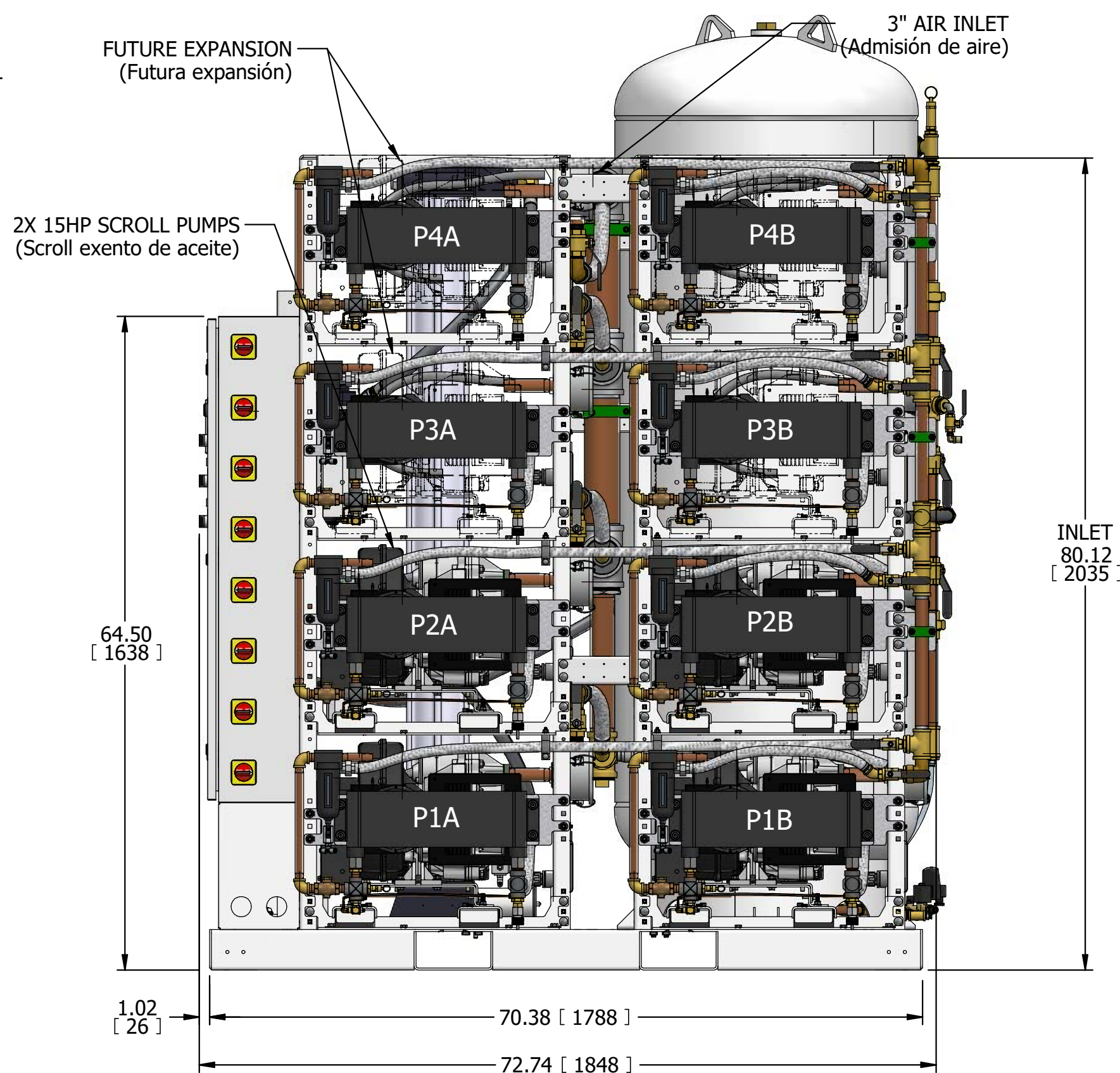
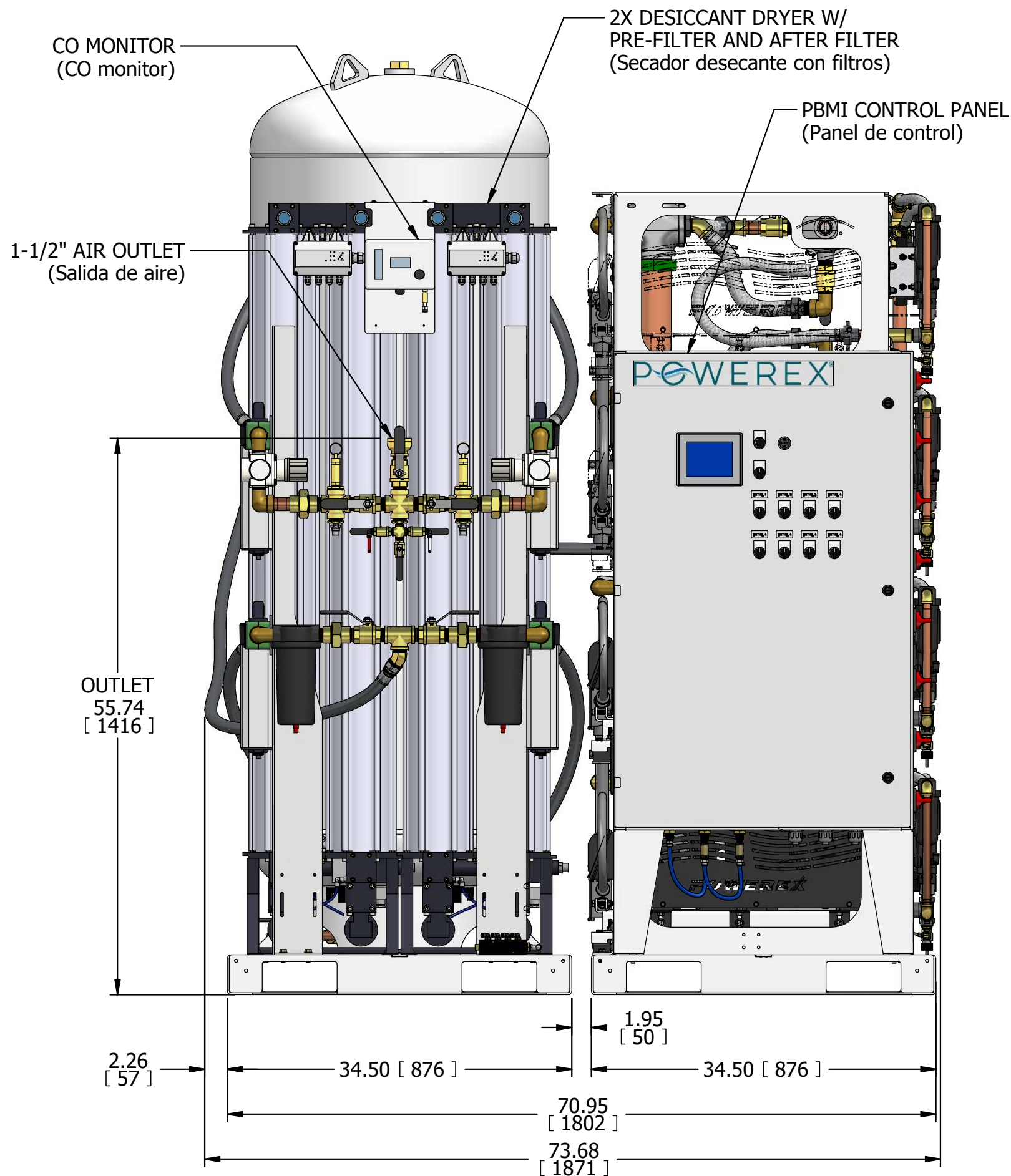
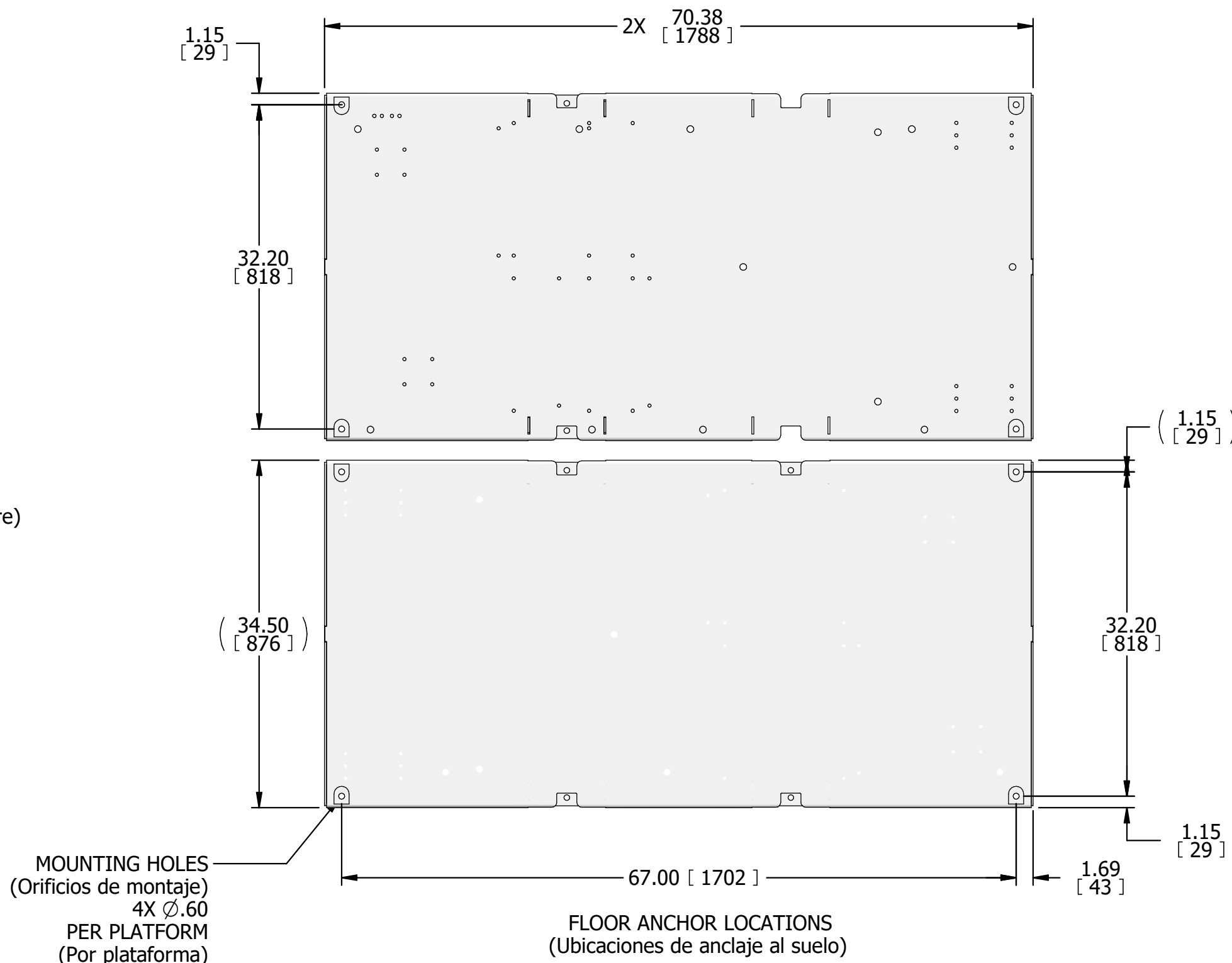



MSD15C6-EX4xP5

15HP MED OPEN SCROLL DUPLEX 240GAL EX4 DUPLEX DESICCANT

REVISION NOTES	
DESCRIPTION	ECN
REVISED CONTROL PANEL FACE	PXEC0318

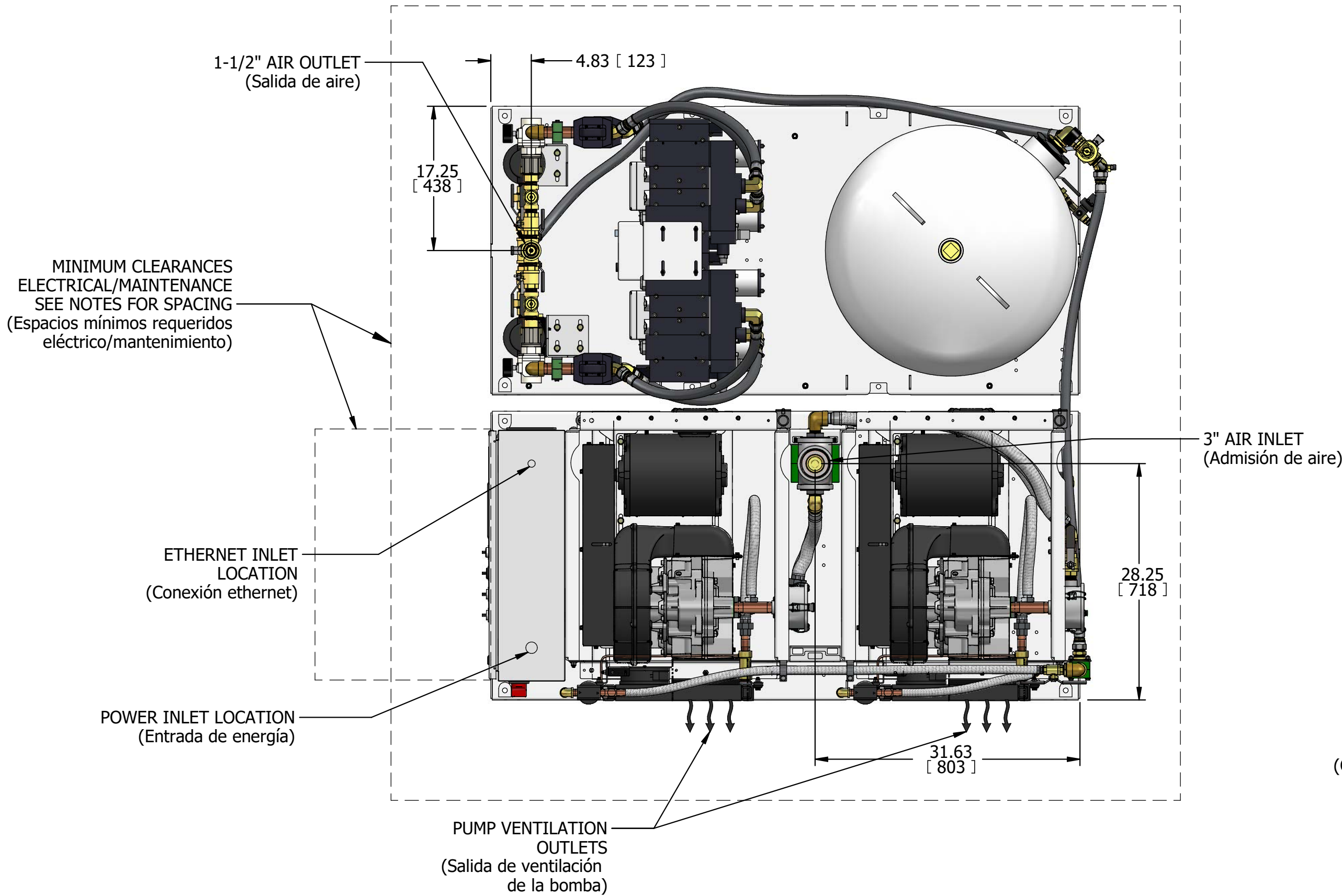

- NOTES:
1. APPROXIMATE WEIGHT: 3806 LBS
 2. A MINIMUM OF 24" CLEARANCE IS RECOMMENDED AROUND THE ENTIRE UNIT. CLEARANCE FOR CONTROL PANEL SHALL BE IN ACCORDANCE WITH OSHA/NEC REQUIREMENTS.
 3. SYSTEM DESIGNED IN COMPLIANCE WITH U.S. NFPA '99, LATEST EDITION.
 4. INLET AND OUTLET FLEX HOSES ADD 12" TO CONNECTIONS
1. PESO APROXIMADO: 1727 kg
2. RECOMENDAMOS QUE EL ESPACIO LIBRE ALREDEDOR DEL SISTEMA DEBE SER COMO MÍNIMO 24" (600 mm). EL ESPACIO LIBRE PARA EL TABLERO DE CONTROL DEBERÁ ESTAR DE ACUERDO CON LOS REQUERIMIENTOS DE SEGURIDAD DE LA OSHA/NEC.
3. EL SISTEMA ESTÁ DISEÑADO DE ACUERDO CON LA ÚLTIMA EDICIÓN DE LA NFPA 99.
4. LAS MANGUERAS FLEXIBLES DE ENTRADA Y SALIDA AÑADEN 12" A LAS CONEXIONES

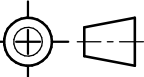
SYMBOL LEGEND				DRAWINGS TO CONFORM TO CURRENT REVISION OF THE FOLLOWING AND STANDARDS: Y14.3 MULTI- AND SECTIONAL VIEW DWGS Y14.5 DIMENSIONING AND TOLERANCING	
TOLERANCES, UNLESS OTHERWISE SPECIFIED		SOLID EDGE DRAWING		DO NOT SCALE DRAWING FOR DIMENSIONS	
0.XX ± .05					
0.XXX ± .005		INCH [mm]		SHEET 1 OF 2	
± .1				SCALE 1:1	
CONFIDENTIAL DISCLOSURE This drawing is the property of POWEREX INC. of HARRISON HOLDINGS, INC. and is subject to return on demand. Its contents are confidential and must not be copied or submitted to outside parties for use or examination.		AUTHOR		REVISION DATE	
		Kumar, Sreehari		07/23/25	
ITEM ID: 761706				APPROVALS ON FILE IN PLM	
MATERIAL: ---					
DESCRIPTION 15HP MED OPEN SCROLL DUPLEX 240GAL EX4 DUPLEX DESICCANT				DRAWING NUMBER MSD15C6-EX4xP5	
				REVISION B	
				SIZE D	

POWEREX

MSD15C6–EX4xP5

15HP MED OPEN SCROLL DUPLEX 240GAL EX4 DUPLEX DESICCANT

REVISION NOTES	
DESCRIPTION	ECN
REVISED CONTROL PANEL FACE	PXEC0318

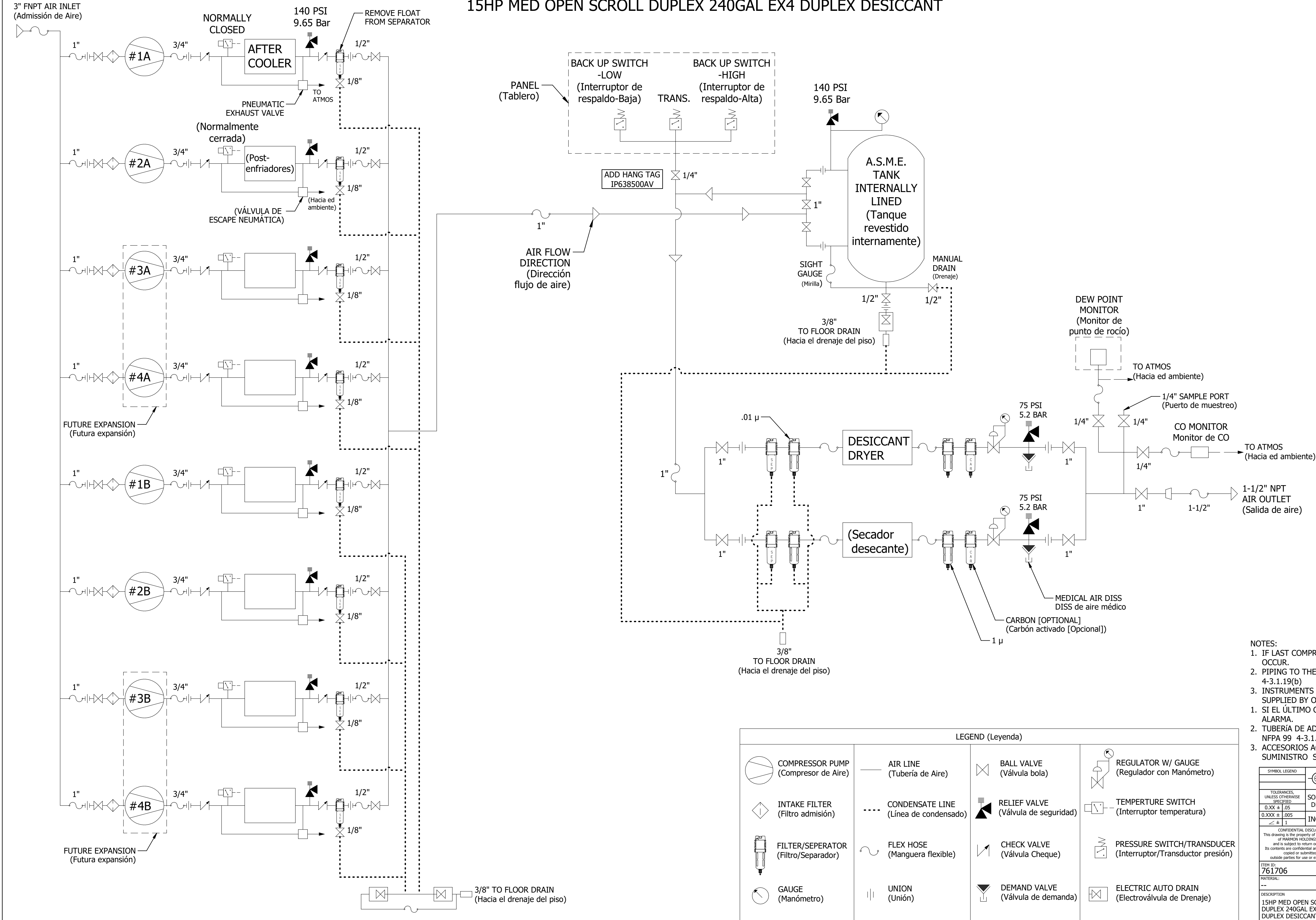

- NOTES:
- IF LAST COMPRESSOR PLEX IS ENERGIZED ALARM WILL OCCUR.
 - PIPING TO THE INLET BY OTHERS AND TO BE PER NFPA99 4-3.1.19(b)
 - INSTRUMENTS DOWN STREAM OF SOURCE VALVE TO BE SUPPLIED BY OTHERS.
 - SI EL ÚLTIMO COMPRESOR SE ACTIVA, SE PRODUCIRÁ UNA ALARMA.
 - TUBERÍA DE ADMISION POR OTROS Y CUMPLIR CON LA NFPA 99 4-3.1.19 (b).
 - ACCESORIOS AGUAS-ABAJO DE LA VALVULA DE CORTE DE SUMINISTRO SERAN INSTALADOS POR OTROS.

LEGEND (Leyenda)			
COMPRESSOR PUMP (Compresor de Aire)	AIR LINE (Tubería de Aire)	BALL VALVE (Válvula bola)	REGULATOR W/ GAUGE (Regulador con Manómetro)
INTAKE FILTER (Filtro admisión)	CONDENSATE LINE (Línea de condensado)	RELIEF VALVE (Válvula de seguridad)	TEMPERTURE SWITCH (Interruptor temperatura)
FILTER/SEPARATOR (Filtro/Separador)	FLEX HOSE (Manguera flexible)	CHECK VALVE (Válvula Cheque)	PRESSURE SWITCH/TRANSDUCER (Interruptor/Transductor presión)
GAUGE (Manómetro)	UNION (Unión)	DEMAND VALVE (Válvula de demanda)	ELECTRIC AUTO DRAIN (Electroválvula de Drenaje)

SYMBOL LEGEND		DRAWINGS TO CONFORM TO CURRENT REVISION OF THE FOLLOWING AND STANDARDS:	
		Y14.3 MULTI AND SECTIONAL VIEW DWGS	
		Y14.5 DIMENSIONING AND TOLERANCING	
TOLERANCES, UNLESS OTHERWISE SPECIFIED	SOLID EDGE DRAWING	DO NOT SCALE DRAWING FOR DIMENSIONS	
0.XX ± .05		SHEET 2 OF 2	SCALE 1:1
0.XXX ± .005	INCH [mm]		
± .1			
CONFIDENTIAL DISCLOSURE This drawing is the property of POWEREX INC. of HARRISON HOLDINGS, INC. and is subject to return on demand. Its contents are confidential and must not be copied or submitted to outside parties for use or examination.		AUTHOR	REVISION DATE
DRAWING NUMBER		REVISION	SIZE
DESCRIPTION			
15HP MED OPEN SCROLL DUPLEX 240GAL EX4 DUPLEX DESICCANT			
DRAWING NUMBER		REVISION	SIZE
MSD15C6–EX4xP5		B	D
MATERIAL:		APPROVALS ON FILE IN PLM	
761706		Kumar,Sreehari 07/23/25	
150 PRODUCTION DRIVE, HARRISON OH 45030 (888) 766-7979		POWEREX®	